

# Data Preprocessing

---

Workshop on RNA-Seq

**Roy Francis and Nima Rafati**

NBIS, SciLifeLab

# Raw data

- Raw count table

| ##                     | DSSd00_1 | DSSd00_2 | DSSd00_3 | DSSd07_1 | DSSd07_2 | DSSd07_3 |
|------------------------|----------|----------|----------|----------|----------|----------|
| ## ENSMUSG000000102693 | 0        | 0        | 0        | 0        | 0        | 0        |
| ## ENSMUSG000000064842 | 0        | 0        | 0        | 0        | 0        | 0        |
| ## ENSMUSG000000051951 | 0        | 1        | 2        | 0        | 3        | 2        |
| ## ENSMUSG000000102851 | 0        | 0        | 0        | 0        | 0        | 0        |
| ## ENSMUSG000000103377 | 0        | 0        | 0        | 0        | 0        | 0        |
| ## ENSMUSG000000104017 | 0        | 0        | 0        | 0        | 0        | 0        |

- Metadata

| ##          | SampleName | SampleID     | No | Model | Day | Group | Replicate |
|-------------|------------|--------------|----|-------|-----|-------|-----------|
| ## DSSd00_1 | DSSd00_1   | KI_PC1606_01 | 1  | DSS   | 0   | day00 | 1         |
| ## DSSd00_2 | DSSd00_2   | KI_PC1606_02 | 2  | DSS   | 0   | day00 | 2         |
| ## DSSd00_3 | DSSd00_3   | KI_PC1606_03 | 3  | DSS   | 0   | day00 | 3         |
| ## DSSd07_1 | DSSd07_1   | KI_PC1606_13 | 13 | DSS   | 7   | day07 | 1         |
| ## DSSd07_2 | DSSd07_2   | KI_PC1606_14 | 14 | DSS   | 7   | day07 | 2         |
| ## DSSd07_3 | DSSd07_3   | KI_PC1606_15 | 15 | DSS   | 7   | day07 | 3         |


# Filtering

- Remove genes and samples with low counts

```
cf1 <- cr[rowSums(cr>0) >= 3, ] # Keep rows/genes that have at least one read
cf2 <- cr[rowSums(cr>2) >= 3, ] # Keep rows/genes that have at least three reads
cf3 <- cr[rowSums(edgeR::cpm(cr)>5) >= 3, ] # need at least three samples to have cpm >
```

count/read per million (cpm/rpm): a normalized value for sequencing depth.

- Inspect distribution

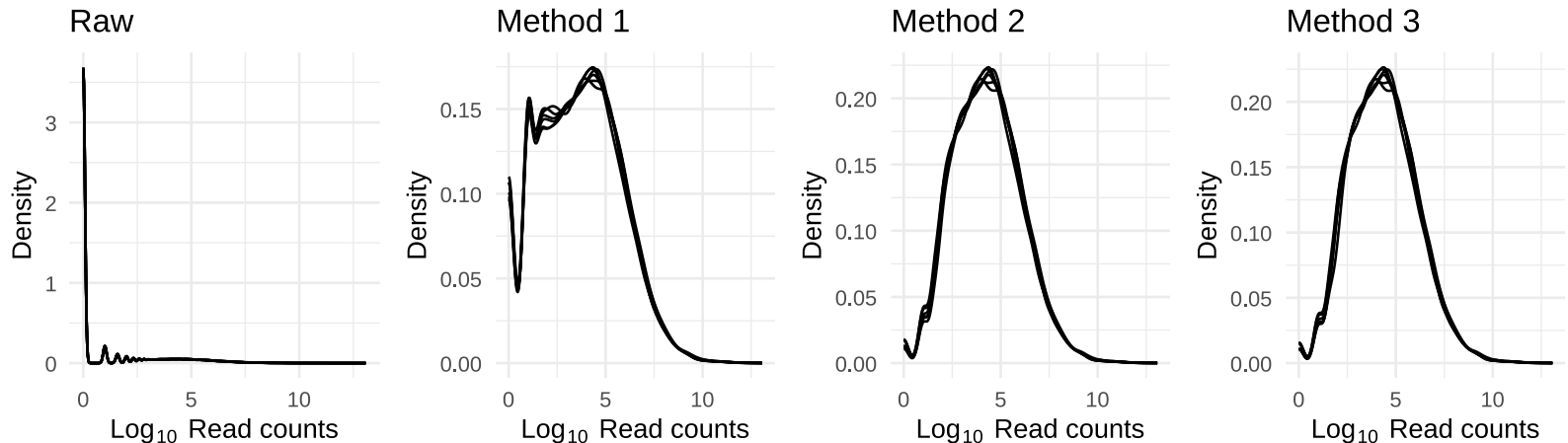


- Inspect the number of rows (genes) available after filtering

```
## Raw: 55487, Method 1: 16099, Method 2: 12656, Method 3: 12496
```

# Normalisation

- Removing technical biases in sequencing data (e.g. sequencing depth and gene length)
- Make counts comparable across samples

- Control for compositional bias

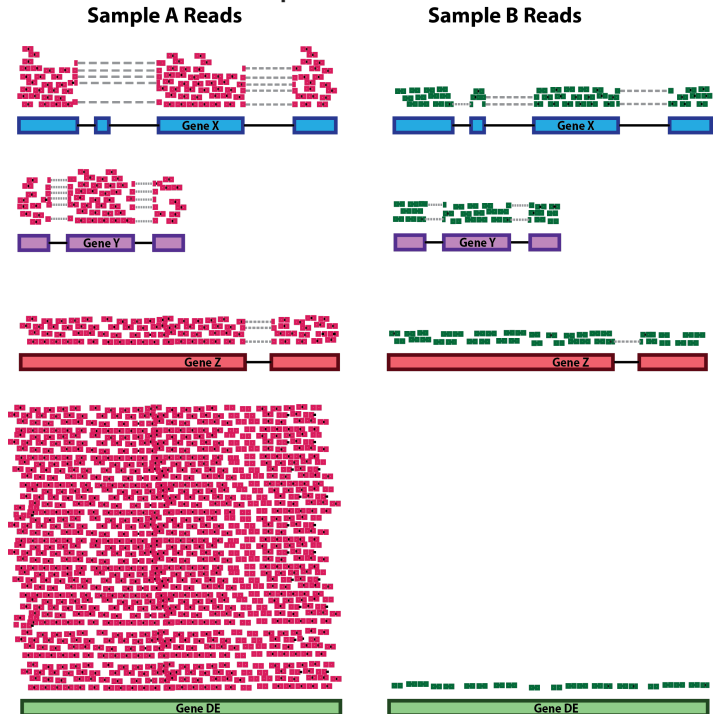

```
##      A B A_tc B_tc
## x  20 6 0.33 0.38
## y  25 6 0.42 0.38
## z  15 4 0.25 0.25
```

```
##      A B A_tc B_tc
## x   20 6 0.12 0.33
## y   25 6 0.16 0.33
## z   15 4 0.09 0.22
## de 100 2 0.62 0.11
```

# Normalisation

- Make counts comparable across features (genes). It can be useful for gene to gene comparisons.

## Sample A Reads



| ##   | counts | gene_length | norm_counts |
|------|--------|-------------|-------------|
| ## x | 50     | 10          | 5           |
| ## y | 25     | 5           | 5           |

- Bring counts to a human-friendly scale

# Normalisation

## Normalisation by library size

- Assumes total expression is the same under different experimental conditions
- Methods include TC, RPKM, FPKM, TPM
- RPKM, FPKM and TPM control for sequencing depth and gene length
- Total number of RPKM/FPKM normalized counts for each sample will be different, therefore, you cannot compare the normalized counts for each gene equally between samples.
- TPM is proportional to RPKM and enables better comparison between samples because total per sample sums to equal value

| ##     | A  | B  | len  | A_rpm   | B_rpm   | A_rpk    | B_rpk    | A_rpk | B_rpk | A_tpm   | B_tpm   |
|--------|----|----|------|---------|---------|----------|----------|-------|-------|---------|---------|
| ## x   | 20 | 6  | 2000 | 408163  | 222222  | 204081.5 | 111111.0 | 10.00 | 3.0   | 493827  | 153846  |
| ## y   | 25 | 6  | 4000 | 510204  | 222222  | 127551.0 | 55555.5  | 6.25  | 1.5   | 308642  | 76923   |
| ## z   | 4  | 15 | 1000 | 81633   | 555556  | 81633.0  | 555556.0 | 4.00  | 15.0  | 197531  | 769231  |
| ## sum | 49 | 27 | 7000 | 1000000 | 1000000 | 413265.5 | 722222.5 | 20.25 | 19.5  | 1000000 | 1000000 |

rpm = cpm.

## Normalisation by distribution

- Assumes technical effects are same for DE and non-DE genes
- Assumes number of over and under-expressed genes are roughly same across conditions
- Corrects for compositional bias
- Methods include Q, UQ, M, RLE, TMM, MRN
- `edgeR::calcNormFactors()` implements TMM, TMMwsp, RLE & UQ
- `DESeq2::estimateSizeFactors()` implements relative log expression (RLE)
- Does not correct for gene length
- `geTMM` is gene length corrected TMM

```
##      A  B  len  ref A_ratio B_ratio      A_mrn      B_mrn
## x 20  6 2000 10.95    1.83    0.55 10.928962 10.90909
## y 25  6 4000 12.25    2.04    0.49 13.661202 10.90909
## z  4 15 1000  7.75    0.52    1.94  2.185792 27.27273
```

## Normalisation by testing

- A more robust version of normalisation by distribution
- A set of non-DE genes are detected through hypothesis testing
- Tolerates a larger difference in number of over and under expressed genes between conditions
- Methods include PoissonSeq, DEGES

## Normalisation using Controls

- Assumes controls are not affected by experimental condition and technical effects are similar to all other genes
- Useful in conditions with global shift in expression
- Controls could be house-keeping genes or spike-ins
- Methods include RUV, CLS

## Stabilizing variance

- Variance is stabilised across the range of mean values
- Methods include VST, RLOG, VROOM
- For use in exploratory analyses. Not for DE.
- `vst()` and `rlog()` functions from *DESeq2*
- `voom()` function from *Limma* converts data to normal distribution



## Recommendations

- Most tools use a mix of many different normalisations
- For DGE using DGE R packages (DESeq2, edgeR, Limma etc), use **raw counts**
- For visualisation (PCA, clustering, heatmaps etc), use VST or RLOG
- For own analysis with gene length correction, use TPM (maybe geTMM?)
- Custom solutions: spike-ins/house-keeping genes



# Thank you. Questions?

R version 4.1.3 (2022-03-10)

Platform: x86\_64-pc-linux-gnu (64-bit)

OS: Ubuntu 18.04.6 LTS

Built on : 🏠 13-Mar-2023 at 🕒 20:05:19

2023 • SciLifeLab • NBIS